




Ergo Table - Height Adjustable Desks with Tilttable Tabletop

Product details

 Table type	B: Whole tabletop tilts
 Operation	Handle
 Dimensions	H: 56-90 cm x W: 90 cm x D: 60 cm (Inner width: 73.6 cm)



Description

The Ergo Table has a tiltable tabletop and manual height adjustment. You can adjust the table height with a handle (56-90 cm), and since the table has no cross bars, wheelchair users can easily get close access to the table.

The tiltable tabletop (0-54°) makes it easy to reach the tabletop, and it also shortens the visual distance to the documents on table. The adjustable tabletop has a matte surface, therefore no annoying glare.

If required, you can use arm supports to relieve your forearms and provide a slight body support for a more balanced sitting position.

Specifications


Item number	20-70011
Power	Manual
Height adjustment	34 cm (56-90 cm)
Tilttable tabletop	Type A: Without tilt. Type B-D: 0-54 ° tilttable tabletop
Operation	Handle
Max. Load on the table (static)	100 kg

Item No: 20-70011

Ergo Table - Height Adjustable Desks with Tilttable Tabletop

Max. Load during lifting (dynamic)	50 kg
Max. load on tilted tabletop (static)	20 kg (type B-D)
Max. load capacity on ruler	8 kg
Materials	Steel St. 37. Various plastic components. Tabletop: Laminate
Surface treatment	Powder coating. Chromate.
Standard colour	Light grey
CE certified	

Line drawing

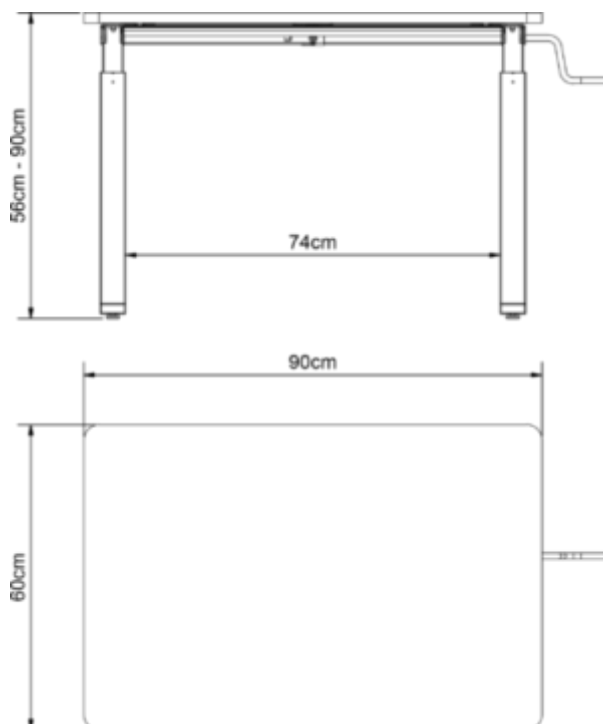

Accessories

Ergo Table - Height Adjustable Desks with Tilttable Tabletop




Item no:

Arm Support

Extra under arm support - can easily be attached to the table and adjusted, including bag and hook. Max user load 25 kg per arm.


Choose variant

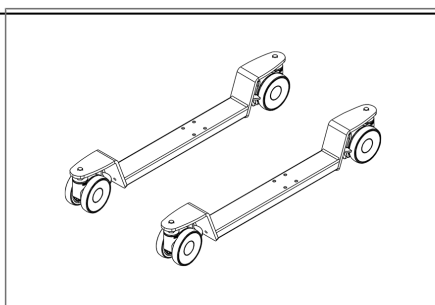

Item no: 20-70502

Factory-installed Wheels - Ergo


Factory-installed brake wheels (set) that make it easy to move the table around. Does not increase the table height. Ordered with the table.

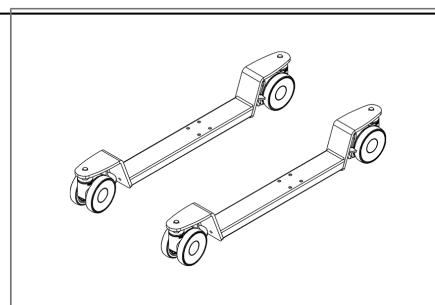

Item no: 20-70516

Wheels For Reordering - Ergo


Brake wheels for retrofitting (set) that make it easy to move the table around. Does not increase the table height.

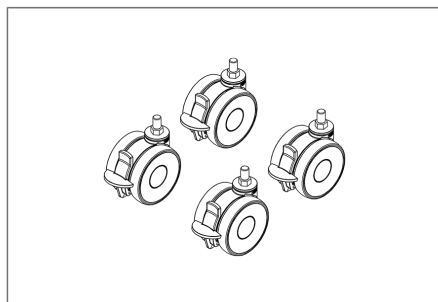

Item no: 20-50340

Brake Wheels With Height Increase

Brake wheels (4 pcs.) that increase the table height by +8 cm.

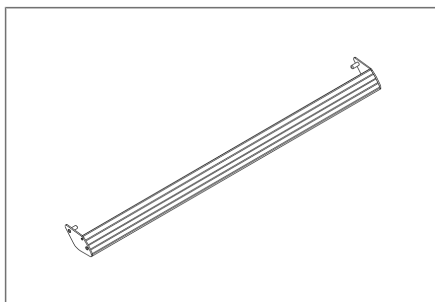


Item no:

Supporting Ruler

Support ruler that can be placed across the entire tabletop to hold paper or a book. Available in different lengths depending on the table model and length of the tabletop. Can only be used on tabletops with an inclined position.

Choose variant

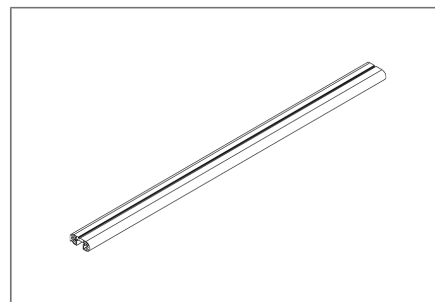

Item no:

Rail For Arm Support

Rail is mounted under the tabletop on the long side, allowing the use of Arm support (1 piece). Used only for table model A = Tabletop without tilt.

Choose variant